


Using EZCalour

EZCalour is a point and click GUI for the Calour microbiome analysis package.

EZCalour can be used to read, process, and plot interactive heatmaps from microbiome experiments.

General concepts

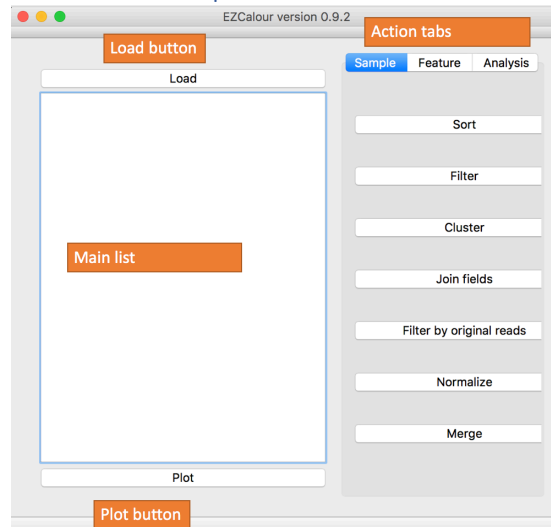

Each dataset is called an experiment. All experiments are displayed in the **main list** in the EZCalour window. Following each processing step, a new experiment is created. Right clicking on an experiment enables deleting from memory and saving of the experiments.

On the right-hand side, there are three **action tabs** - for processing **samples**, **features** and **analysis**. Commands from a given tab relate to the appropriate axis (i.e. Filter from the sample tab filters samples, whereas Filter from the features tab filters features). Commands work on the selected experiment from the main list.


In order to plot the interactive heatmap of an experiment, it needs to be selected from the main list, and then press the **"Plot"** button (located at the lower left side)

Loading data

EZCalour works with microbiome BIOM tables, metabolomics MZMine2 tables, or any CSV text file.

Besides the main table, ezcalour can also load a tab-separated mapping file, containing information about each sample.

In order to load an experiment, click on the "Load" button (located at the top left side).

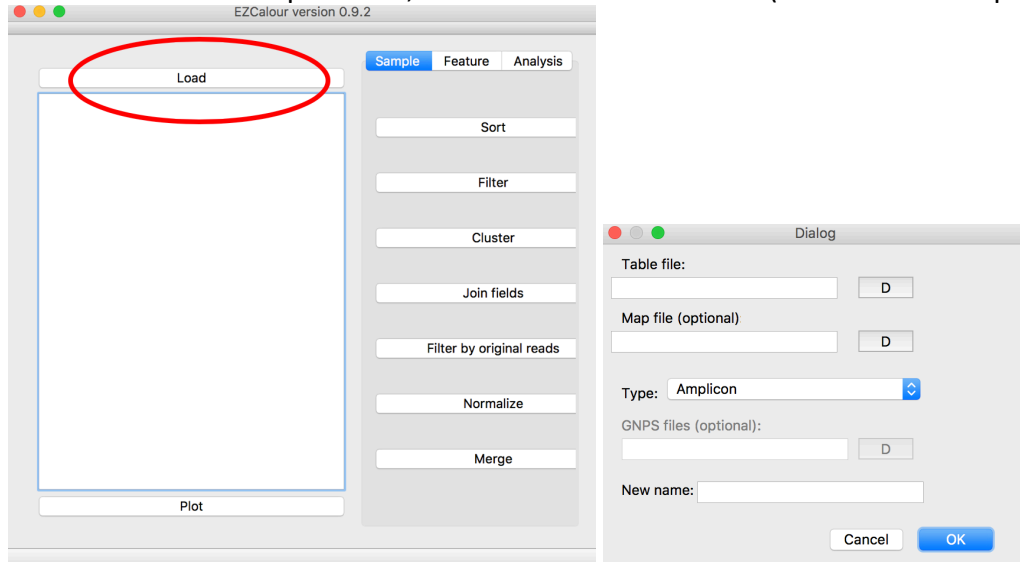

Mandatory fields:

- "Table file" : name of the biom or mzmine2 table (can click the "D" button for GUI file selection)
- "Type": the type of the table file:
 - "Amplicon" for a microbiome biom table. When loading, the table is normalized by TSS to 10000 reads/sample. Samples with <1000 reads are dropped.
 - "MZMine2" for an MZMine2 metabolomics table
 - "TSV" for a general tab separated table (Each sample is a column, each feature is a row)

Optional fields:

- "Map file" : name of the sample TSV mapping file
- "GNPS file" : For mass-spec, the per-metabolite info file (see [here](#))
- "New name" : the name for the experiment in the main list (defaults to the table file name)

Processing data


Sample tab

Buttons from this tab affect the samples of the experiment. Button actions are performed on the selected experiment from the main list. Each action generates a new experiment in the main list.

Sort

Sort the samples according to the selected field.

Sorting is conservative, meaning samples with same value in the field retain the previous order. In order to sort by two fields (i.e. "Disease" and "Day" within each disease), sort first by the second field (i.e. "Day") and then by the first (i.e. "Disease").

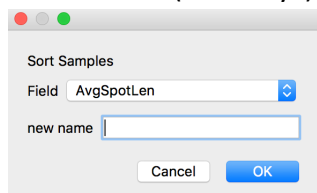

Mandatory fields:


- "Field": select the sample metadata field to sort by

Optional fields:

- "New name" : the name for the experiment in the main list (defaults to the table file name)

Filter

Keep or remove samples with specified mapping file field values

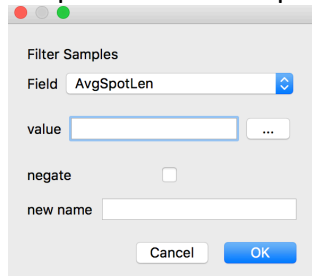

Mandatory fields:

- "Field": select the sample metadata field to sort by
- "Value": the values to filter for the field.

Optional fields:


- "neagte": if checked, remove samples with the selected values, otherwise keep samples with selected values
- "New name" : the name for the experiment in the main list (defaults to the table file name)

Cluster

Cluster the samples by putting similar samples next to each other

Join Fields

Create a new metadata field by joining the values of two fields.

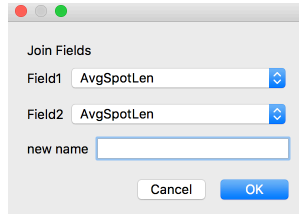A dialog box titled "Join Fields" with a light gray background. It contains two dropdown menus, "Field1" and "Field2", both set to "AvgSpotLen". Below them is a text input field labeled "new name" which is empty. At the bottom are "Cancel" and "OK" buttons.

Mandatory fields:

- "Field1", "Field2": The two fields to join (new field will be field1-field2).


Optional fields:

- "New name" : the name for the experiment in the main list.

Filter by original reads

Throw away samples with < threshold original reads (before normalization).

Used to get rid of samples with a small number of reads

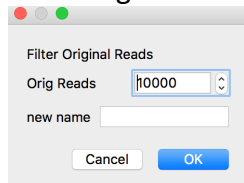A dialog box titled "Filter Original Reads" with a light gray background. It contains a spinner control labeled "Orig Reads" set to "10000". Below it is a text input field labeled "new name" which is empty. At the bottom are "Cancel" and "OK" buttons.

Mandatory fields:

- "Orig. reads": the minimal number of reads in the sample in order to keep

Optional fields:


- "New name" : the name for the experiment in the main list.

Normalize

TSS (total sum scaling) normalization of the reads per sample.

When loading a microbiome table, samples are automatically scaled to 10000 reads/sample.

NOTE: This is not rarefaction. So features can have fractional reads.

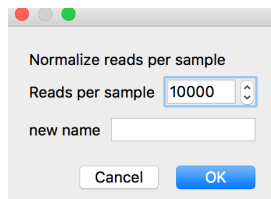A dialog box titled "Normalize reads per sample" with a light gray background. It contains a spinner control labeled "Reads per sample" set to "10000". Below it is a text input field labeled "new name" which is empty. At the bottom are "Cancel" and "OK" buttons.

Mandatory fields:

- "Reads per sample": The sum of reads per sample to normalize to


Optional fields:

- "New name" : the name for the experiment in the main list.

Merge

merge samples having the same value in the selected field.

Samples can be merged using mean, sum or randomly choosing one of the samples with the value.

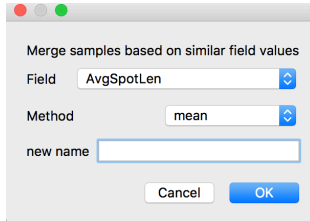

Mandatory fields:

- "Field": The field to merge samples sharing the same value
- "Method": how to merge the samples with same value. Options are:
 - "mean": new sample contains the mean of each of the features from all the samples with the value
 - "random": new sample contains a randomly chosen sample with the value
 - "sum": new sample contains the sum of each of the features from all the samples with the value

Optional fields:

- "New name" : the name for the experiment in the main list.

Feature tab


Buttons from this tab affect the features of the experiment. Button actions are performed on the selected experiment from the main list. Each action generates a new experiment in the main list.

Cluster

Cluster the features by putting similar behaving features next to each other.

Each feature is normalized to mean=0, std=1 and then Euclidian distance is used for clustering, so similar behaving features over the samples will be close to each other, without dependence on the absolute level of the features.

Can also remove low abundance features (to speed up the clustering).

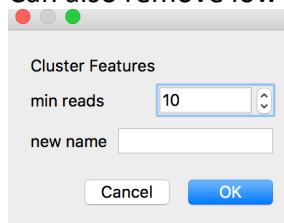A dialog box titled "Cluster Features" with a title bar containing red, yellow, and green window control buttons. It contains two input fields: "min reads" with a numeric value of 10 and a spin button, and "new name" with an empty text box. At the bottom are "Cancel" and "OK" buttons.

Mandatory fields:


- "min reads": Keep only features with total reads over all samples \geq min reads.

Optional fields:

- "New name" : the name for the experiment in the main list.

Filter min reads

Keep only "interesting features" that have enough total reads (over all samples)

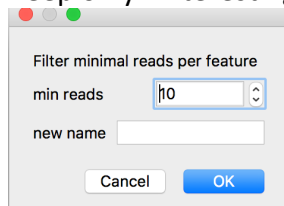A dialog box titled "Filter minimal reads per feature" with a title bar containing red, yellow, and green window control buttons. It contains two input fields: "min reads" with a numeric value of 10 and a spin button, and "new name" with an empty text box. At the bottom are "Cancel" and "OK" buttons.

Mandatory fields:

- "min reads": Keep only features with total reads over all samples \geq min reads.


Optional fields:

- "New name" : the name for the experiment in the main list.

Filter taxonomy

Keep (or remove) features matching a given taxonomy string.

NOTE: This requires the biom table to contain taxonomy information.

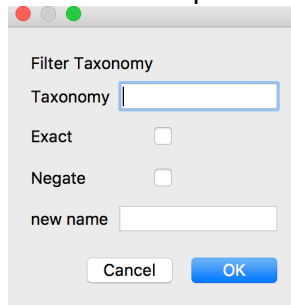A dialog box titled "Filter Taxonomy" with a light gray background. It contains a text input field labeled "Taxonomy" with a blue border. Below it are two checkboxes: "Exact" and "Negate", both of which are unchecked. At the bottom is a text input field labeled "new name". At the very bottom are two buttons: "Cancel" and "OK".

Mandatory fields:


- "Taxonomy": A partial or complete taxonomy string. Needs to match the taxonomy embedded in the biom table.

Optional fields:

- "Exact": Check to filter only taxonomy strings fully matching the string. Uncheck to allow partial matches (i.e. "roteo" will match "p__Proteobacteria")
- "Negate": Check to remove matching features, uncheck to keep matching features.
- "New name" : the name for the experiment in the main list.

Filter fasta

Keep (or remove) features appearing in a fasta file

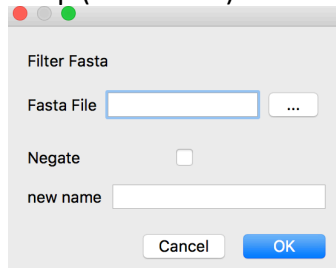A dialog box titled "Filter Fasta" with a light gray background. It contains a text input field labeled "Fasta File" with a blue border, followed by a button with three dots "...". Below it is a checkbox labeled "Negate" which is unchecked. At the bottom is a text input field labeled "new name". At the very bottom are two buttons: "Cancel" and "OK".

Mandatory fields:


- "Fasta file": Name of the fasta file containing the feature IDs (usually sequences or qiime2 hashes).

Optional fields:

- "Negate": Check to remove matching features, uncheck to keep matching features.
- "New name" : the name for the experiment in the main list.

Filter prevalence

Keep features appearing at least in a given fraction of the samples (common features)

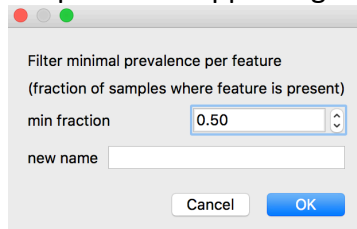

A dialog box titled "Filter minimal prevalence per feature" with the subtitle "(fraction of samples where feature is present)". It contains a "min fraction" input field with a value of 0.50 and a "new name" text field. At the bottom are "Cancel" and "OK" buttons.

Mandatory fields:


- "min fraction": the minimal fraction of the samples where the feature is present.

Optional fields:

- "New name" : the name for the experiment in the main list.

Filter mean

Keep features with a large enough mean frequency (over all samples)

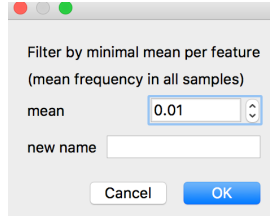

A dialog box titled "Filter by minimal mean per feature" with the subtitle "(mean frequency in all samples)". It contains a "mean" input field with a value of 0.01 and a "new name" text field. At the bottom are "Cancel" and "OK" buttons.

Mandatory fields:


- "mean": the minimal mean frequency (over all samples) for features to be kept.

Optional fields:

- "New name" : the name for the experiment in the main list.

Sort abundance

Order the features based on their mean frequency over (an optional subset of) samples.

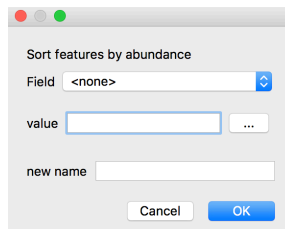

A dialog box titled "Sort features by abundance". It contains a "Field" dropdown menu currently set to "<none>", a "value" input field, and a "new name" text field. At the bottom are "Cancel" and "OK" buttons.


Optional fields:

- "Field": <none> to sort based on frequency over all samples, otherwise sort based on frequency only based on samples matching "value" in "Field".
- "value": The value to use for the sample subset.
- "New name" : the name for the experiment in the main list.

Collapse taxonomy

Merge features based on their taxonomy (summing the frequencies of all features with the same taxonomy).

NOTE: taxonomy information must be embedded in the biom table.

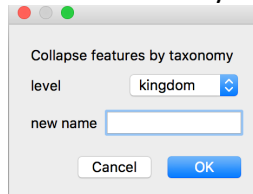

Collapse features by taxonomy

level kingdom

new name

Cancel OK

Mandatory fields:

- "level": The taxonomic level to merge by (i.e. Phyla, Genus, etc.)

Optional fields:

- "New name" : the name for the experiment in the main list.

Analysis tab


Buttons from this tab are used for statistical analysis of the experiment. Button actions are performed on the selected experiment from the main list. Most actions generate a new experiment in the main list.

Diff. abundance

Find statistically significant features separating two subsets of samples from the experiment. All the tests are non-parametric with dsFDR control.

Result is a new experiment, containing only the significant features, sorted by the effect size (topmost is the highest in group1, bottommost is the highest in group2). In order to view the significant bacteria, just plot this new experiment.

NOTE: There is no special separation marker between the last feature higher in group1 and the first feature higher in group2, since we believe if you don't see it by eye, it's not that interesting.

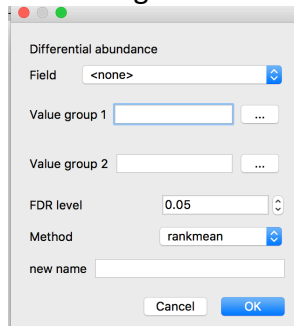

Mandatory fields:

- "Field": Name of the field used to separate the samples into two groups.
- "Value group 1": the value used for samples in group1. Can select multiple values using the "... " button.
- "Value group 2": the value used for samples in group2. Can select multiple values using the "... " button. If no values are selected, all samples not in group1 are used for group2.


Optional fields:

- "FDR level": The upper bound on features where the null hypothesis holds (similar to Benjamini-Hochberg q value).
- "Method": The statistic used for the comparison between the two groups. Options are:
 - "rankmean": rank each feature over the samples and compare the mean between the two groups. Robust to outliers.
 - "mean": just compare the means of the two groups for each feature. Can be strongly affected by outliers.
 - "binary": Compare the fraction of samples where the feature is present.
- "New name" : the name for the experiment in the main list.

Correlation

Find statistically significant features correlated with a metadata field. All the tests are non-parametric with dsFDR control.

Result is a new experiment, containing only the significant features, sorted by the effect size (topmost is the highest correlated feature, bottommost is the highest anti-correlated). In order to view the significant bacteria, just plot this new experiment.

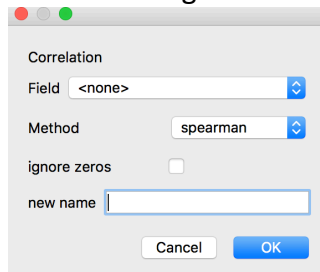

Mandatory fields:

- "Field": Name of the field used to look for correlation.

Optional fields:


- "Method": The statistic used for the comparison between the two groups. Options are:
 - "spearman": spearman rank correlation (rank samples and each feature over samples and then calculate pearson correlation).
 - "pearson": pearson correlation between feature frequency and metadata field.
 - "ignore zeros": For each feature, calculate the correlation only on samples where the feature is present (NOTE: much slower).
- "New name" : the name for the experiment in the main list.

Enrichment

Find dbBact terms enriched in features from the diff. abundance test.

Shows the dbBact terms which are significantly more associated with features higher in group2 or group2 (from the diff. abundance test).

NOTE: can only be used on experiments which are the result of diff. abundance.

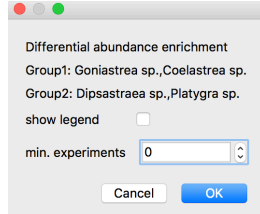

Differential abundance enrichment

Group1: Goniastrea sp., Coelastrea sp.

Group2: Dipsastraea sp., Platygra sp.

show legend ☐

min. experiments

Cancel OK

Optional fields:

- “show legend”: check this box to show the name of each group in the resulting bar plot
- “min. experiments”: Keep only terms appearing in at least min. experiments different experiments (i.e. if set to 2, ignore terms coming only from one experiment).

Result is a bar plot showing terms enriched in either of the two groups.

Example:

